


GABA Research ELISA

Featuring Patent Pending Extraction Technology



Kit Details

Cat No.	Sample Size	Standards	Sensitivity	Format
BA E-2500	300 µl serum, or any biological fluid	0 / 15 - 1 500 ng/ml	7.5 ng/ml	 96 Wells

Features

LDN has developed an ELISA for the quantitative determination of GABA in different biological samples. The determination of GABA has the following features:

- Easy sample extraction step featuring a patent pending technology that avoids the use of ion-exchange columns
- No tedious derivatization
- Freedom from organic solvents
- Highly specific for GABA (no cross reactions with Glutamate or Glutamine)
- Broad dynamic measuring range
- Freedom from interferences

New GABA ELISA

based on an unique high throughput extraction procedure

Faassen, H.J.R. van, and Manz B., Labor Diagnostika Nord, Nordhorn, Germany

Abstract

Introduction: The new GABA (γ -aminobutyric acid) ELISA is especially designed to measure GABA in different biological fluids (eg. brain and spinal cord homogenates, cell culture). With this new assay improved sensitivity is achieved and less sample volume can be used. To achieve all this a new extraction procedure was developed. Instead of using the laborious ion exchange column procedure, the ion exchange resin was coated to a microtiter plate (patent pending). It showed that after the coating, the ion exchange kept its unique characteristics. The advantage of using an extraction plate instead of a column: no "dead" volume, faster throughput, higher sensitivity. For the subsequent ELISA procedure a highly specific GABA antibody is used. A proof of principle was performed with rat brain and spinal cord tissue homogenates. Samples of rat brain and spinal cord were homogenized and spiked with known concentrations of GABA and recovery and linearity of the GABA ELISA was tested.

Results: GABA standards were treated like samples and corresponded to 250 - 25000 pg/well. Only 10 μ L homogenate was needed for the analysis. Recovery from brain homogenate (cortex) was 83 - 100 %, for spinal cord homogenate (lumbar region) 85 - 92 %. Linearity for both samples was 98 - 115 %. Different brain and spinal cord regions were analyzed (brain: cortex and hippocampus, spinal cord: cervical, thoracic and lumbar region). Results were in line with previous HPLC results on the same kind of samples.

Conclusion: A sensitive GABA ELISA with a new extraction procedure was developed. Based on this straightforward method high throughput testing of samples is possible. The results show that it is possible to measure GABA in different kind of homogenates.

Extraction

Sample Preparation (Flexible Volumes)

1. GABA is extracted by ion exchange chromatography. Microtiter plates (48 Well) are coated with a cation exchange resin (patent pending).
2. After washing away unbound molecules, the bound GABA is released. Whole extraction procedure takes 30 minutes only.
3. In a simple derivatization step GABA is coupled to a protein. This is needed for the ELISA

Samples


- Brain homogenates
- Spinal cord homogenates
- Cell culture supernatants

Immunoassay

1. 50 μ L derivatized sample
2. 50 μ L specific antiserum
3. Incubation overnight at 4 - 8 $^{\circ}$ C
4. Three wash cycles with 300 μ L wash buffer
5. 100 μ L Enzyme conjugate
6. Incubation for 30min at RT (shaker)
7. Three wash cycles with 300 μ L wash buffer
8. 100 μ L Substrate
9. Incubation for 20-30min at RT (shaker)
10. 100 μ L Stop solution
11. Measure at 450nm

Validation

Typical Standard Curve

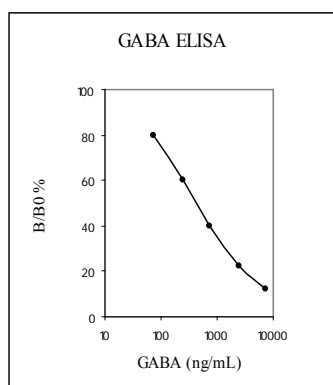

Recoveries

Sample	Recovery
Motor cortex	Mean recovery: 88 % SD=7.9; CV=9.0 R ² =0.999
Spinal cord Lumbar region	Mean recovery: 87 % SD=3.6; CV=4.1 R ² =0.997
Sample	Linearity
Hippocampus	Mean linearity: 105 % SD=8.2; CV=7.8 R ² =0.997

Measured Rat Samples*

Tissue sample	Measured conc. (pmol/mg wet weight tissue)
Motor cortex	1137 \pm 65
Somatosensory cortex	1379 \pm 115
Hippocampus	1719 \pm 219
Cervical sc region	484 \pm 65
Thoracic sc region	395 \pm 19
Lumbar I & II sc region	540 \pm 80
Lumbar III - VI sc region	737 \pm 32
Sacral sc region	687 \pm 48

*Samples were kindly provided by E. Ziemińska, Nencki Institute, Poland

World-Wide:

Labor Diagnostika Nord GmbH & Co KG Am Eichenhain 1, 48531 Nordhorn, Germany Tel: +49 5921 8197-101
Fax: +49 5921 8197-222 Web: www.ldn.de e-mail: info@ldn.de

North America:

Rocky Mountain Diagnostics, Inc. 703 N Tejon Suite C - Colorado Springs, CO 80903 Tel: 1-877-477-0039
Fax: 1-719-477-0307 Web: www.rmdiag.com e-mail: info@rmdiag.com